
# QUICK START GUIDE

US EN

For a full explanation of all features and instructions, please refer to the User's Guide. (available for download from [www.motorolastore.com](http://www.motorolastore.com)).

## 1. Overview of the control panel

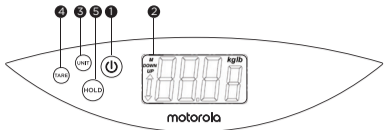- 1. Power ON/OFF button**
- 2. LCD Display**
- 3. UNIT button** Select the unit of measurement: kilograms (kg) or pounds (lb).
- 4. TARE button** Zeros the weight of a towel placed on the platform.

### 5. HOLD unit

The function HOLD allows to stabilize the weight on the screen when the baby is quite still and his weight has not changed, i.e. the hands of an adult while they weigh the child.

Press HOLD while the baby is on the platform alone and still, so the machine stabilises the detected weight (it will flash on the screen). When the display stops flashing, the weight of the baby will be displayed alternating with the up/down difference from the last weight.

## 2. Install batteries

The battery operates with 4 AAA alkaline batteries type LR03 (1.5V), not included.

- The battery housing is on the bottom of the weighing scale, carefully pry open the battery cover on the underside of the Scale.
- Install four (4) new "AAA" batteries (LR03) into the battery compartment, making sure to match the positive and negative ends of the batteries to the correct contacts as indicated in the battery compartment.
- Replace the battery cover after insert the batteries.
- Press the Power ON/OFF button to turn on the weighing scale.

### Note

- Batteries are to be installed with the correct polarity.
- Do not mix batteries of different types or mix old and new batteries, replace all 4 batteries when empty with new batteries.
- If the scale is to be stored unused for a long period of time, the batteries should be removed.
- The battery contacts are not to be short-circuited.

## 3. Using the baby scale

### 3.1 Turn On/Off the baby scale

Press the Power ON/OFF button to turn On/Off the weighing scale.

### 3.2 Select Weighing Unit

Press the UNIT button on the control panel to select unit of measurement, kilograms (kg) or pounds (lb).

### 3.3 Tare Weight

The weighing scale zeros automatically when turned on and displays the digit "0.00", if you weight the baby with a soft cloth or towel, it need to tare weight.

- To add a towel or cloth when the machine is already turned on.
- Press the TARE button to zero the weight of soft cloth or towel placed on the platform, and weight that the digit "0.00" is displayed.

### 3.4 Weight Tracker

The digital baby scale allows to follow the evolution of the baby's weight during his growth, with a maximum load of 30 kg and a division of 10 g.

#### 3.4.1 First weighing

- Press Power ON/OFF button to turn on the baby scale, when the digits "0.00" are displayed, the machine is ready for weighing (If necessary, use the function TARE as previously described section 4.3).
- Place the baby on the weighing scale and the weight will be displayed on the screen.
- Press the HOLD button to stabilize the readings , and the weight data will be saved in the memory of baby scale.

#### 3.4.2 Following weighings

- To turn on the baby scale by pressing Power ON/OFF button.
- Place the baby on the weighing scale when the digit "0.00" is displayed.
- The weight of the baby will appear on the screen, press the HOLD button to stabilize the readings.
- The measured weights stop flashing and following weight data will display as alternating.
  - Weight of baby
  - Weight Change (Compare with last weighing)

- if it is weight loss, the symbol and weight change are display on the screen (Fig.1).

- if it is weight gain, the symbol and weight change are display on the screen (Fig.2).

- No weight change was detected (Fig.3).

- The weight data will be saved in the memory of baby scale, the weighing scale with turn off after 20 seconds if no any operation.

### NOTE

The weight data saved in scale memory will be erased when batteries are removed or replaced, the weighing scale will start again from "First Weighing" to record the weight data.

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